


ERGONOMIC POSITIONING:
DUAL AUTO-REGRESSION

The dual regression system of the **GALAXY PLUS PEDIATRIC** bed allows the backrest section to be adjusted (11cm), reducing bedsores to a minimum.

IMPROVE THE PATIENTS' POSITION
AND RECOVERY TIME

The cardiac chair position is fundamental for patients suffering from heart and respiratory problems, aiding their recovery. Medical personnel may sit the patient in our new bed at the touch of a button, using the automated positioning system of the **GALAXY PLUS PEDIATRIC** beds.

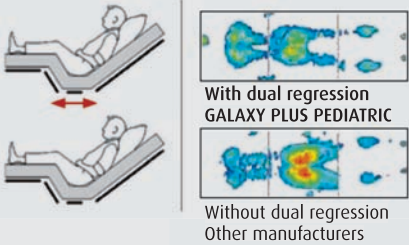

With dual regression
GALAXY PLUS PEDIATRIC

Without dual regression
Other manufacturers



POSITIONING AND COMFORT



Satellite control panel
Clear pictograms for patients and medical staff.



Nurses control panel
A panel for the sole use of medical personnel that allows different bed functions to be activated or overridden. It also incorporates certain preD programmed functions at the touch of a button (Trendelenburg, CPR, cardiac seat and memory position)



Patient control panel
The patient's control panel includes two functions that allow the bed to be positioned in the cardiac seat position and fully horizontal.

EASE OF CLEANING

Our high quality products reduce cleaning time by up to 30%. Smooth surfaces and removable panels allow easy and fast cleaning, and disinfection.




BRAKE PEDAL ALARM

When the bed is connected to the electricity supply, this alarm alerts medical staff when the brake system is activated. The bilateral height-adjustment pedal allows the bed height to be adjusted to the most adequate position in order to reduce the risk of injury. As an added safety function, the pedal is automatically blocked when not in use.



EASY MOVEMENT

The fifth wheel facilitates the rotation of the bed, reducing the effort when moving it.

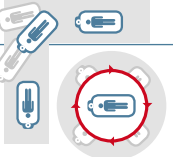

GALAXY PLUS PEDIATRIC Bed

Quality support for the youngest

BASIC EQUIPMENT

1. Completely motorised bed with scissor system.
2. Maximum permitted load: 230 Kg.
3. Bed base made up of 4 independent ABS modules.
4. Integrated bed extension (+15cm).
5. Patient control panel.
6. Nurse's control panel.
7. GO safety button.
8. Autoregression (11 cm backrest / 7 cm pelvis).
9. Long-lasting rechargeable battery.
10. 150 mm diameter castors with fairing.
11. Protective bumpers.
12. Lower fairing in ABS.
13. Four supports for serum posts and lifting pole.
14. Cable support.
15. Mattress retainers.
16. Central braking system.
17. Steering castor.
18. Back rest angle indicator and TR/ATR.
19. Bilateral supports for accessories.
20. 8 hooks for restraint straps.

ACCESSORIES

1. Holder for urine bag (REF. CHA20).
2. Lifting pole (REF. CHA17).
3. IV Drip pole (REF CHA15).

OPTIONAL EQUIPMENT

1. Metal fold-down side rails (REF. CHV50).
2. Independent ABS side rails (REF. CHV53).
3. Side railings with controls (REF CHV84).
4. Bilateral height adjustment pedal with blocking function (REF. CHV83).
5. Fifth wheel (REF. CHV 24).
6. Brake alarm activated / deactivated (REF. CHA30) .
7. Adjustable lamp (REF. CHA28).
8. Satellite control panel (REF. CHV81).
9. X-Ray cassette holder (REF. CHV26).

ELECTRICAL SPECIFICATIONS

Power supply	230/110V 50/60Hz
Maximum current intake	1.6A, 370 Volt.
Liquid ingress protection indicator	IP 66
Protection class	Class I
Electrical shock protection	Type B

TECHNICAL DATA

External dimensions	190 x 86 cm.
Patient surface	177x70 cm.
Height adjustment	37,5 - 77,5 cm.
Backrest tilting	70°
Legrest tilting	45°
Backrest autoregression plane	11 cm.
Trendelemburg / antitrendelemburg	+ - 16°
Weight without accessories	150 kg.
Maximum working load	230 kg.
Bed extension	+ - 15 cm.



BUREAU VERITAS



ECO-LABEL
ISO 14001



IBV



CE

www.medical-iberica.com

Lanzahita, 6 - Pol. Ind. Albarreja - 28946 Fuenlabrada, Madrid (Spain) - T. 91 208 94 00 - F. 91 208 94 01

NOTE: The dimensions and characteristics of the products displayed may be subject to change without notice.

medisa
cares for you

Clinical Furniture
GALAXY PLUS PEDIATRIC Bed

GALAXY PLUS PEDIATRIC Bed

Quality support for the youngest

The GALAXY PLUS PEDIATRIC has been developed to meet the highest demands of a standard hospital environment.



- It reduces the physical effort to move patients
- It helps reduce bedsores
- Innovative concept of side handrails that extend along the length of the bed
- With the new GALAXY PLUS PEDIATRIC bed system, we help speed up the patients' recovery thanks to an advanced ergonomic design that increases their autonomy and mobility.

NEW SAFETY REQUIREMENTS

According to standard EN 60 601-2-52 for electric beds, handrails should minimize the risk of trapping the patient's head, neck or limbs to avoid possible fatal accidents.

Ergonomic head and footboard panels help manoeuvre the bed. The materials used meet all the necessary requirements for hygiene and maintenance.

Bed base made up of 4 independent HPL modules

TOTAL PROTECT side rails: offering perfect protection against falls. Compliant with standard EN 60 601-2-52



Trendelenburg /Antitrendelenburg: +16°

Patient Control panel

CPR position control which allows backrest to be positioned horizontally

Integrated bed extension

Protective bumpers at each of the corners

Universal support bar for accessories

Centralised braking system and power steering facility